

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# NFKB p65 (phospho-Ser536) polyclonal antibody

Catalog: BS66162 Host: Rabbit Reactivity: Human, Mouse, Rat

#### **BackGround:**

NF-kappa-B is a pleiotropic transcription factor present in almost all cell types and is the endpoint of a series of signal transduction events that are initiated by a vast array of stimuli related to many biological processes such as inflammation, immunity, differentiation, cell growth, tumorigenesis and apoptosis. NF-kappa-B is a homo- or heterodimeric complex formed by the Rel-like domain-containing proteins RELA/p65, RELB. NFKB1/p105, NFKB1/p50, REL and NFKB2/p52 and the heterodimeric p65-p50 complex appears to be most abundant one. The dimers bind at kappa-B sites in the DNA of their target genes and the individual dimers have distinct preferences for different kappa-B sites that they can bind with distinguishable affinity and specificity. Different dimer combinations act as transcriptional activators or repressors, respectively. NF-kappa-B is controlled by various mechanisms of post-translational modification and subcellular compartmentalization as well as by interactions with other cofactors or corepressors. NF-kappa-B complexes are held in the cytoplasm in an inactive state complexed with members of the NF-kappa-B inhibitor (I-kappa-B) family. In a conventional activation pathway, I-kappa-B is phosphorylated by I-kappa-B kinases (IKKs) in response to different activators, subsequently degraded thus liberating the active NF-kappa-B complex which translocates to the nucleus. NF-kappa-B heterodimeric p65-p50 and p65-c-Rel complexes are transcriptional activators. The NF-kappa-B p65-p65 complex appears to be involved in invasin-mediated activation of IL-8 expression. The inhibitory effect of I-kappa-B upon NF-kappa-B the cytoplasm is exerted primarily through the interaction with p65. p65 shows a weak DNA-binding site which could contribute directly to DNA binding in the NF-kappa-B complex. Associates with chromatin at the NF-kappa-B promoter region via association with DDX1.

#### **Product:**

Rabbit IgG, 1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

### **Molecular Weight:**

~ 70 kDa

#### **Swiss-Prot:**

Q04206

#### **Purification&Purity:**

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

#### **Applications:**

WB 1:500 - 1:1000

IHC 1:50 - 1:200

IF/ICC 1:50 - 1:200


#### Storage&Stability:

Store at  $4 \, \mathbb{C}$  short term. Aliquot and store at  $-20 \, \mathbb{C}$  long term. Avoid freeze-thaw cycles.

#### **Specificity:**

NFKB p65 (phospho-Ser536) polyclonal antibody detects endogenous levels of NFKB p65 protein only when phosphorylated at Ser536.

#### **DATA:**

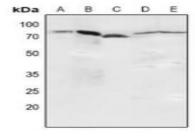

Western blot analysis of NF-kappaB p65 (Phospho-S536) expression in Hela (A), HEK293T (B), NIH3T3 (C), mouse kidney (D), H9C2 (E) whole cell lysates.

#### Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841

# Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,


P. R. China.

Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Tel: 0086-025-68037686
Fax: 0086-025-68035151




# PRODUCT DATA SHEET

Bioworld Technology,Inc.

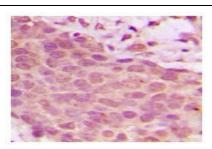

Immunohistochemical analysis of NF-kappaB p65 (Phospho-S536) staining in human breast cancer formalin fixed paraffin embedded tissue section.

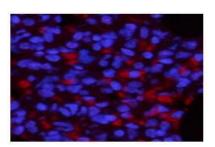

 $Immun of luorescent\ analysis\ of\ NF-kappa B\ p65\ (Phospho-S536)\ staining\ in\ rat\ lung.$ 

#### Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151