

MBOAT4 polyclonal antibody

Catalog: BS31641

Host: Rabbit

Reactivity: Human, Mouse

BackGround:

Membrane-bound ghrelin O-acyltransferase MBOAT4: Catalyzes ghrelin acylation at 'Ser-3' using preferentially octanoyl-CoA, hexanoyl-CoA and decanoyl-CoA as acyl-CoA donors leading to ghrelin activity. In vitro also uses acyl-CoA donors of different lengths from short-chain (C2) to long-chain fatty acids (C16) knowing that acyl-CoA donors from butanoyl-CoA (C4) to dodecanoyl-CoA (C12) are more efficient compared to longer acyl-CoA donors, such as myristoyl-CoA (C14) and palmitoyl-CoA (C16) that are not efficient

Product:

1mg/ml in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% sodium azide.

Molecular Weight:

~49.7 kDa

Swiss-Prot:

Q96T53

Purification&Purity:

The antibody was purified by immunogen affinity chromatography.

Applications:

WB (1/500 - 1/1000), IHC (1/50 - 1/200)



Storage&Stability:

Store at 4 °C short term. Aliquot and store at -20 °C long term. Avoid freeze-thaw cycles.

Specificity:

Recognizes endogenous levels of MBOAT4 protein.

DATA:

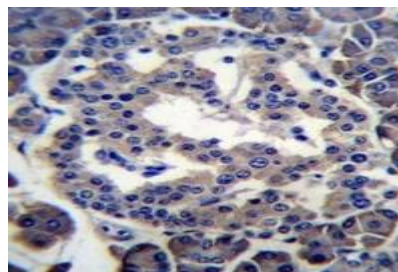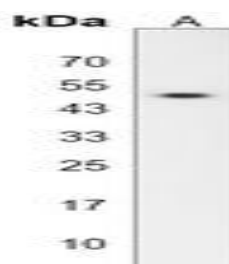

Immunohistochemical analysis of OCT2 staining in human Hodgkin's lymphoma formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

Immunohistochemical analysis of OCT2 staining in human tonsil formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park, MN 55416, USA.

Email: info@bioworld.com

Tel: 6123263284

Fax: 6122933841

Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046, P. R. China.

Email: info@biogot.com

Tel: 0086-025-68037686

Fax: 0086-025-68035151